

EYEPIECE ADAPTERS

PM-ADF

PM-ADP

PM-ADG


OM
SYSTEM

Summary

- These are designed to mount an eyepiece on a microscope for photomicrography with the OM System cameras.

Main Characteristics

- Three types are available in accordance with the type of Olympus eyepieces FK, P, and G in use.
- The Adapter PM-ADF is designed for Photo Eyepiece FK, the Adapter PM-ADP for Viewing Eyepiece P, and the Adapter PM-ADG for Stereoscopic Eyepiece G.

Notes • Related Units

- The letters "FK", "P" or "G" engraved on each eyepiece indicate the eyepiece type.
- Photo eyepieces FK are preferable for photomicrography.
- If picture magnifications (or total magnifications) are the same, optimum structural details of the object can be photographed by means of a higher-objective-and-lower-eyepiece combination.

OLYMPUS

OM SYSTEM EYEPIECE ADAPTERS PM-ADF/PM-ADP/PM-ADG

Eyeiece Adapter PM-ADF




● **Camera Assembly Mount**
Accepts the OM-Mount Photomicro Adapter L, Automatic Exposure Body PM-PBA or Manual Exposure Body PM-PBM.

● **Clamping Screw**


● **Microscope Mount**
Fits into the eyepiece sleeve of a microscope and tightens with clamping screw.

Eyeiece Adapter PM-ADP

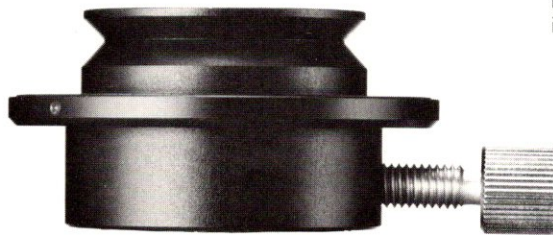


Eyeiece Adapter PM-ADG


Consists of an upper and a lower part. When used with an Eyepiece G10X, insert the G10X into the lower part and screw the upper part on it before mounting the entire unit to the eyepiece sleeve.

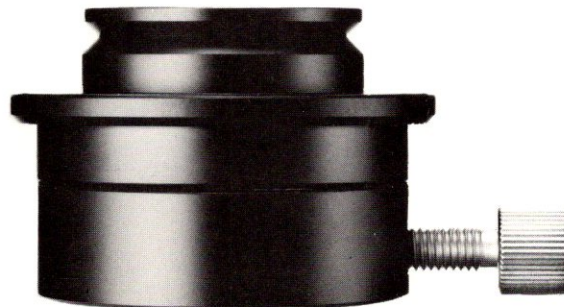

Eyeiece Adapter PM-ADF

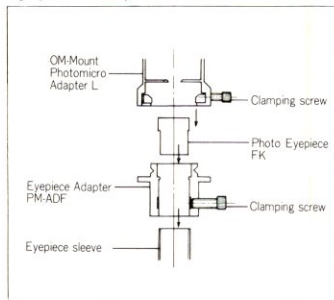

Eyeiece Adapter PM-ADP

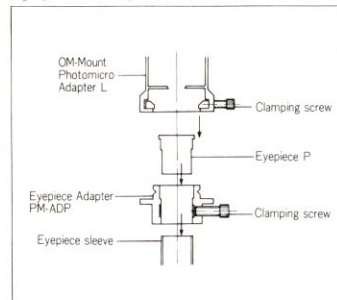

Eyeiece Adapter PM-ADG

